

Visian ICL™ TORIC

eDFU-0042/1

P h a k i c T o r i c I n t r a o c u l a r L e n s

TABLE OF CONTENTS

| | | | |
|---|---|---|---|
| PRODUCT INFORMATION | 1 | SUMMARY OF SAFETY AND CLINICAL PERFORMANCE..... | 2 |
| DEVICE DESCRIPTION..... | 1 | CLINICAL TRIAL TICM..... | 2 |
| INDICATIONS..... | 1 | ADVERSE EVENT REPORTING..... | 3 |
| MODE OF ACTION..... | 1 | HOW SUPPLIED..... | 3 |
| CONTRAINDICATIONS | 1 | EXPIRATION DATE..... | 3 |
| COMPLICATIONS AND ADVERSE REACTIONS..... | 1 | RETURN POLICY FOR STAAR VISIAN TICL..... | 3 |
| PRECAUTIONS..... | 2 | WARRANTY AND LIMITATION OF LIABILITY | 3 |
| CALCULATION OF LENS POWER AND SIZING | 2 | STORAGE | 3 |
| LENS PREPARATION..... | 2 | CAUTION..... | 3 |
| ADMINISTRATION AND INSTRUCTION FOR USE..... | 2 | REFERENCES / BIBLIOGRAPHY | 4 |
| WARNINGS..... | 2 | SYMBOL GLOSSARY | 5 |

DIRECTIONS FOR USE

PRODUCT INFORMATION

Please review this product information completely before performing your initial clinical procedure. All physicians must complete the STAAR Surgical Visian® TICL Physician Certification Program, special attention is placed on sizing methodologies for determination of the Visian ICL overall diameter. Improper ICL size may lead to adverse events ranging from mild to severe.

DEVICE DESCRIPTION


The Visian® Toric Implantable Collamer® Lens (Visian TICL) features a single piece lens design with a central, concave/convex optical zone of 4.7 to 5.8 mm diameter (according to model and diopter). The lens is manufactured in four overall diameters, 11.6, 12.1, 12.6, 13.2 mm to accommodate different eye sizes. The lenses are capable of being folded and implanted through an incision of 3.5 mm or less.

The lenses are manufactured from a proprietary ultraviolet (UV) radiation absorbing polymer containing hydroxyethylmethacrylate (HEMA) and porcine collagen. The 10% UV cut-offs for STAAR's phakic IOL lens family are:

- 377 nm for the thinnest central thickness lens, -5.5 D, and
- 388 nm for the thickest central thickness lens, +10.0 D.

Table 1: VTICH Models

| Brand Name | Model Name | Dioptic Power (D) | Cylindrical power | Overall Diameter (mm) | Optic Diameter (mm) | Haptic Design |
|-------------|------------|-------------------|-------------------|-----------------------|---------------------|---------------|
| Visian TICL | VTICH 11.6 | 0.0 to +10.0 | +0.5 to +6.0 | 11.6 | 4.7 to 5.8 | Flat, plate |
| Visian TICL | VTICH 12.1 | 0.0 to +10.0 | +0.5 to +6.0 | 12.1 | 4.7 to 5.8 | Flat, plate |
| Visian TICL | VTICH 12.6 | 0.0 to +10.0 | +0.5 to +6.0 | 12.6 | 4.7 to 5.8 | Flat, plate |
| Visian TICL | VTICH 13.2 | 0.0 to +10.0 | +0.5 to +6.0 | 13.2 | 4.7 to 5.8 | Flat, plate |

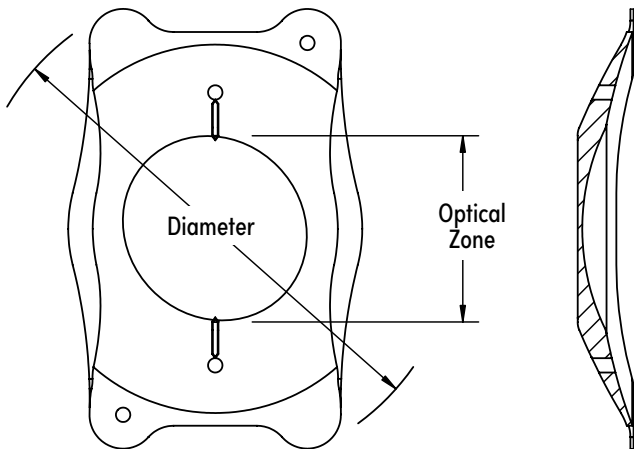

VTICH Diagram

INDICATIONS

Visian® Toric Implantable Collamer® Lenses (Visian TICL) are indicated for use in phakic eye treatment in patients 21-45 years of age for:

- The correction/reduction of hyperopia in patients ranging from +0.5D to +16.0D with less than or equal to 6.0D of astigmatism at the spectacle plane.
- With an anterior chamber depth (ACD) equal to or greater than 3.0 mm, as measured from the corneal endothelium to the anterior lens capsule.

MODE OF ACTION

Visian TICL is intended to be placed entirely within the posterior chamber directly behind the iris and in front of the anterior surface of the natural lens of a phakic patient. When correctly positioned, the lens functions as a refractive element to optically correct/reduce hyperopia with astigmatism.

CONTRAINDICATIONS

Visian TICL is contraindicated in the presence of any of the following circumstances and/or conditions:

1. Patients with low/abnormal corneal endothelial cell density, Fuchs' dystrophy or other corneal pathology
2. Ocular hypertension in either eye
3. Any cataract in the operative eye or non-traumatic cataract in the fellow eye.
4. Persons under the age of 21 years.
5. Primary Open Angle or Narrow Angle Glaucoma.
6. Narrow anterior chamber angles (i.e. less than Grade III as determined by gonioscopic exam).
7. Pregnant or nursing.
8. Previous or pre-existing ocular disease that would preclude post-operative visual acuity of 0.477 logMAR (20/60 Snellen) or better.
9. Patients who are amblyopic or blind in the fellow eye.
10. Implantation of a lens in an eye with an anterior chamber depth (ACD), as measured from the corneal endothelium to the anterior lens capsule, less than 3.0 mm.

COMPLICATIONS AND ADVERSE REACTIONS

Adverse reactions and complications due to, or following surgery and implantation of any Visian TICL may include, but are not limited to: Hyphema, Non-reactive Pupil, Pupillary Block, Additional YAG Iridotomy, Secondary Glaucoma, Cataract, Intraocular Infection, Uveitis/ Iritis, Retinal Detachment, Vitritis, Corneal Edema, Macular Edema, Corneal Decompensation, Over/Under Correction, Significant Glare and/or Halos (under night driving conditions), Hypopyon, Increased Astigmatism, Loss of BSCVA, Rotation/Decentration/Subluxation, IOP Elevation from Baseline, Corneal Endothelial Cell Loss, Iris Pigment Dispersion, Secondary Surgical Intervention to Remove/ Replace/Reposition the Lens, Peripheral Anterior Synechia (PAS), Iris Synechia to Implant, Conjunctival Irritation, Vitreous Loss.

PRECAUTIONS

1. Do not attempt to sterilize.
2. Do not autoclave.
3. The lens should not be exposed to any solutions other than the normally used intraocular irrigating solutions (e.g. isotonic saline, Balanced Salt Solution (BSS), viscoelastic, etc.).
4. The lens should be handled carefully. No attempts should be made to reshape or cut any portion of the lens or to apply undue pressure to the lens optical portion with a sharp object.
5. Do not allow the lens to dry in air. The lens should be stored in sterile BSS solution during surgery.
6. The long-term effect of the lens has not been determined. Therefore, physicians should continue to monitor implant patients postoperatively on a regular basis.
7. Safety and effectiveness of the lens has not been established in patients with: unstable refractive error in either eye, keratoconus, history of clinical signs of iritis/uveitis, synechia, pigment dispersion syndrome, pseudoexfoliation, insulin-dependent diabetes or diabetic retinopathy, history of previous ocular surgery including refractive corneal surgery.
8. Implantation of a lens may result in a decrease in corneal endothelial cell density.

CALCULATION OF LENS POWER AND SIZING

The lens power and size calculation should be performed by the surgeon using the STAAR OCOS Calculation Software. Using the software potentially prevents calculation errors that may result in secondary surgery due to refractive surprise, excessive vaulting, lens rotation, IOP elevation from baseline, etc. During the U.S. FDA trial for the ICM/TICM, the white to white and ACD (from the corneal endothelium to the anterior lens capsule) were used to determine the ICL overall diameter. There are some reports suggesting that white to white corneal measurements do not correlate with sulcus to sulcus measurements. Recent publications indicate that new imaging technologies may provide optimal visualization and measurement of the intraocular dimensions involved in phakic intraocular lens implantation.

LENS PREPARATION

Verify that the level of the liquid fills at least 2/3 of the vial. The thermoformed tray and vial should be opened in a sterile field. Record serial number on operative report to retain traceability of the lens. Remove the aluminum cap and stopper from the vial. Remove the lens from the vial. The lens should not be exposed to a dry environment (air) for more than one minute.

CAUTION: Do not use if package has been opened or damaged.

CAUTION: Do not allow the lens to dry after removal from the glass vial.

ADMINISTRATION AND INSTRUCTION FOR USE

Implantation of the Visian TICL should only be attempted by a surgeon who is highly skilled in the required surgical technique. The following procedure is recommended for implantation of Visian TICL. Two YAG iridotomies, (0.5 to 0.8 mm; placed superiorly, 90 degrees apart) should be performed 1 to 2 weeks prior to surgery with confirmation of patency prior to lens implantation. The patient should be prepared for surgery according to the surgeon's standard operating procedure. Mark the desired axis (TARGET axis) for alignment of Visian TICL. A clear scleral or corneal tunnel wound incision of 3.5 mm or less should be used, followed by filling of the anterior and posterior chamber with an appropriate viscoelastic. The lens is then folded using a MICROSTAAR® injector MSI-PF or MSI-TF with SFC-45 cartridge and injected into the anterior chamber of the eye. Please refer to the MICROSTAAR® injection system Directions for Use for instructions regarding proper loading and injection of the lens using the MICROSTAAR® injection system. Verify correct orientation of the lens in the eye and that the lens is not inverted. If the pupil remains sufficiently dilated, the lens should be well centered and positioned under the iris in front of the natural lens of a phakic patient, so that the footplates are placed in the sulcus. Verify the

Visian TICL alignment marks are in the required axis (TARGET axis) according to the Implantation Orientation Diagram (IOD). Complete removal of the viscoelastic material from the eye must be performed after completion of the surgical procedure and before the eye is closed (without sutures). From this point the operation can proceed according to the surgeon's standard procedure. Dispose of any single use accessories that may have become contaminated with bodily fluids during the procedure as biohazardous waste according to standard surgical biohazard waste disposal procedure. Postoperative medical care of the patient should also follow the surgeon's standard procedure.

WARNINGS

1. Check the label of the lens package for proper lens model and power.
2. Open the package to verify the dioptric power of the lens.
3. Handle the lens by the haptic portion. Do not grasp the optic with forceps and never touch the center of the optic once the lens is placed inside the eye.
4. Complete removal of viscoelastic from the eye after completion of the surgical procedure is essential. STAAR Surgical recommends a low molecular weight 2% hydroxypropyl methylcellulose (HPMC) or dispersive, low viscosity ophthalmic viscosurgical device.
5. STAAR Surgical recommends using the MICROSTAAR® MSI-PF or MSI-TF with SFC-45 cartridge injection systems to insert the lens in the folded state.

NOTE: The primary viscoelastic used during the US FDA clinical trial was a low molecular weight 2% hydroxypropyl methylcellulose preparation.

SUMMARY OF SAFETY AND CLINICAL PERFORMANCE

A Summary of Safety and Clinical Performance (SSCP) for the ICL family of lenses can be found in the European database on medical devices (Eudamed) at <https://ec.europa.eu/tools/eudamed>. The Basic UDI-DI (BUDI-DI) used to search for the for the ICL family of lenses on the website is 764013516ICLGV.

CLINICAL TRIAL TICM

Summary Findings of the Clinical Study:

The TICM Implantable Collamer® Lenses were found to be safe and effective as refractive elements to optically reduce moderate to high myopia with astigmatism.

Table 2: Patient Demographics
210 Eyes Treated of 124 Patients

| Sex | | |
|---------------|-----|-------|
| Male | 55 | 44.4% |
| Female | 69 | 55.6% |
| Ethnic Origin | | |
| Caucasian | 102 | 82.3% |
| Hispanic | 10 | 8.1% |
| Black | 6 | 4.8% |
| Other | 6 | 4.8% |

Mean Age 35 years

Table 3: Adverse events

The adverse events reported in the 210 eyes enrolled during the clinical trial (at any postoperative exam) are presented below

| Adverse Event | N | % |
|---|---|------|
| Raised IOP Requiring Treatment | 1 | 0.5% |
| Pupillary Block | 1 | 0.5% |
| Retinal Detachment | 1 | 0.5% |
| Surgical Reintervention | | |
| TICL Repositioning | 1 | 0.5% |
| Visian TICL Replacement (too long) | 1 | 0.5% |
| Visian TICL Removal (no ICL or IOL replacement) | 3 | 1.4% |
| YAG Iridotomy | 3 | 1.4% |
| BCDVA loss ≥ 2 lines | 3 | 1.5% |
| Anisocoria | 1 | 0.5% |

Table 4: Best Corrected Distance Visual Acuity (BCDVA) with time for patients with preop BCDVA 20/20 or better

| | 6 Months n/N, % | 12 Months n/N, % |
|----------|--------------------|---------------------|
| ≤20/12.5 | 71/155, 45.8% | 72/159, 45.3% |
| ≤20/16 | 141/155, 91.0% | 143/159, 89.9% |
| ≤20/20 | 155/155, 100% | 159/159, 100% |
| ≤20/25 | 155/155, 100% | 159/159, 100% |
| ≤20/40 | 155/155, 100% | 159/159, 100% |

Table 5: Uncorrected Distance Visual Acuity (UCDVA) over time for patients with Preoperative BCDVA 20/20 or Better

| | Preop n/N, % | 6 Months n/N, % | 12 Months n/N, % |
|----------|-----------------|--------------------|---------------------|
| ≤20/12.5 | 0/173, 0% | 41/155, 26.5% | 40/159, 25.2% |
| ≤20/16 | 0/173, 0% | 117/155, 75.5% | 101/159, 63.5% |
| ≤20/20 | 0/173, 0% | 140/155, 90.3% | 142/159, 89.3% |
| ≤20/40 | 0/173, 0% | 155/155, 100% | 159/159, 100% |
| >20/50 | 173/173, 100% | 0/155, 0.0% | 0/159, 0.0% |
| >20/200 | 173/173, 100% | 0/155, 0.0% | 0/159, 0.0% |

Table 6: Manifest Refraction Spherical Equivalent over Time

| | Preop | 1 Week | 1 Month | 3 Months | 6 Months | 12 Months |
|-------------------------------|-----------------|---------------|---------------|---------------|---------------|----------------|
| N (eyes) | 210 | 205 | 200 | 191 | 182 | 194 |
| Mean Spherical Equivalent (D) | -9.38 | 0.02 | 0.13 | 0.13 | 0.11 | 0.03 |
| SD | 2.67 | 0.45 | 0.43 | 0.39 | 0.49 | 0.46 |
| Range (D) | -19.50 to -2.38 | -1.50 to 1.38 | -1.63 to 1.75 | -1.25 to 1.25 | -1.75 to 2.63 | -2.25 to ±1.00 |

Table 7: Manifest Refraction Cylinder over Time

| | Preop | 1 Week | 1 Month | 3 Months | 6 Months | 12 Months |
|-------------------------------|--------------|--------------|--------------|--------------|--------------|--------------|
| N (eyes) | 210 | 205 | 200 | 191 | 182 | 194 |
| Mean Spherical Equivalent (D) | 1.95 | 0.50 | 0.50 | 0.52 | 0.45 | 0.52 |
| SD | 0.84 | 0.54 | 0.49 | 0.49 | 0.45 | 0.48 |
| Range (D) | 1.00 to 4.00 | 0.00 to 3.00 | 0.00 to 3.00 | 0.00 to 3.00 | 0.00 to 2.00 | 0.00 to 3.00 |

Table 8: Rotation of the TICL Between Visits (from direct observation of TICL)

| Rotation | 1 Day – 1 Week n/N, % | 1 Week – 1 Month n/N, % | 1 Month – 3 Months n/N, % | 3 Months – 6 Months n/N, % | 6 Months – 12 Months n/N, % |
|----------|--------------------------|----------------------------|------------------------------|-------------------------------|--------------------------------|
| ≤ 5° | 118/121, 97.5% | 148/155, 95.5% | 141/148, 95.3% | 133/136, 97.8% | 132/140, 94.3% |
| ≤ 10° | 121/121, 100% | 155/155, 100% | 147/148, 99.3% | 135/136, 99.3% | 137/140, 97.9% |

ADVERSE EVENT REPORTING

Adverse Reactions and/or potentially sight-threatening complications that may reasonably be regarded as lens related must be reported immediately to STAAR Surgical and the competent authority of the Member State where the patient is established. This information is being requested from surgeons in order to document potential long-term effects of Visian® TICL implantation

HOW SUPPLIED

Visian TICL is supplied sterile and non pyrogenic in a sealed vial containing BSS. The vial is sealed within a thermoformed tray placed in a box with labels and product information. Sterility is assured until the expiration date indicated on package label, if the tray and vial seal are not punctured or damaged. Visian TICL is steam sterilized. A Patient Implant Card, implant card instructions and labels are supplied in the unit package. This card which includes a link to important safety information regarding the implanted lens, should be completed by the healthcare provider and given to the patient to keep as a permanent record of the implant, and a resource to show to any eye care practitioner seen in the future.

EXPIRATION DATE

The expiration date on the device package is the sterility expiration date. This device must not be used past the indicated sterility expiration date.

RETURN POLICY FOR STAAR VISIAN TICL

Contact STAAR Surgical. Visian TICL must be returned dry. Do not attempt to re-hydrate the lens.

WARRANTY AND LIMITATION OF LIABILITY

STAAR Surgical warrants that reasonable care was taken in making this product. STAAR Surgical shall not be responsible for any incidental or consequential loss, damage, or expense which arises directly or indirectly from the use of this product. To the extent permitted by law, STAAR Surgical's sole liability from any and all causes pursuant to Visian TICL shall be limited to the replacement of Visian TICL which is returned to and found to be defective by STAAR Surgical. This warranty is in lieu of and excludes all other warranties not expressly set forth herein, whether expressed or implied by operation of law or otherwise, including but not limited to, or any implied merchantability or fitness for use.

STORAGE

Store the lens at room / ambient temperature.

CAUTION

- Do not autoclave the lens. Do not store the lens at temperatures greater than 40°C. Do not freeze. If temperature requirements are not met, return the lens to STAAR Surgical.
- STAAR Surgical Visian TICLs and disposable accessories are packaged and sterilized for single use only. Cleaning, reuse and/or resterilization are not applicable to these devices. If one of these devices were reused after cleaning and/or resterilization it is highly probable that it would be contaminated and the contamination could result in infection and/or inflammation.

REFERENCES / BIBLIOGRAPHY:

1. Bylsma SS, Zalta AH, Foley E, Osher RH. Phakic Posterior Chamber Intraocular Lens Pupillary Block. *J Cataract Refract Surg.* 2002;28:2222-2228.
2. Sanders DR, Actual and Theoretical Risks for Visual Loss Following Use of the Implantable Contact Lens for Moderate to High Myopia. *J Cataract Refract Surg.* 2003;29:1323-1332.
3. Sanders DR, Vukich JA, Doney K, Gaston M. U.S. Food and Drug Administration Clinical Trial of the Implantable Contact Lens for Moderate to High Myopia, *Ophthalmology* 2003;110:255-266.
4. Sanders DR: Postoperative Inflammation after Implantation of the Implantable Contact Lens. *Ophthalmology* 2004;110:2335-41.
5. Edelhauser HF, Sanders DR, Azar R, Lamielle H. Corneal endothelial assessment after ICL implantation. *J Cataract Refract Surg* 2004;30:576-583.
6. Sanders DR, Doney K, Pico M. U.S. Food and Drug Administration Clinical Trial of the Implantable Collamer Lens for Moderate to High Myopia, Three-Year Follow-up. *Ophthalmology* 2004;111:1683-1692.
7. Lege BAM, Haigis W, Neuhann TF, Bauer MH. Age-Related Behavior of Posterior Chamber Lenses in Myopic Phakic Eyes During Accommodation Measured by Anterior Segment Partial Coherence Interferometry. *J Cataract Refract Surg.* 2006;32:999-1006.
8. Sanders DR. Matched Population Comparison of the Visian Implantable Collamer Lens and Standard LASIK for Myopia of -3.00 to -7.88 Diopters. *J Refract Surg.* 2007;23:537-553.
9. Sanders DR, Sarver EJ. Standardized Analyses of Correction of Astigmatism With the Visian Toric Phakic Implantable Collamer Lens. *J Refract Surg.* 2007;23:649-660.
10. Schallhorn S, Tanzer D, Sanders DR, Sanders ML. Randomized Prospective Comparison of Visian Toric Implantable Collamer Lens and Conventional Photorefractive Keratectomy for Moderate to High Myopic Astigmatism. *J Refract Surg.* 2007;23:853-867.
11. Sanders DR, Schneider D, Martin R, Brown D, Dulaney D, Vukich J, Slade S, Schallhorn S. Toric Implantable Collamer Lens for Moderate to High Myopic Astigmatism. *Ophthalmology* 2007;114:54-61.
12. Sanders DR, Anterior Sub-capsular Opacities and Cataracts 5 Years After Surgery in the Visian Implantable Collamer Lens FDA Trial. *J Refract Surg* 2008;24:566-570.
13. Sanders DR, Sanders ML. Comparison of the Toric Implantable Collamer Lens and Custom Ablation LASIK for Myopic Astigmatism. *J Refract Surg.* 2008;24:773-778.
14. Sanders DR, Bernitsky DA, Harton PJ, Rivera RP. The Visian Myopic Implantable Collamer Lens Does Not Significantly Affect Axial Length Measurement With the IOL Master. *J Refract Surg.* 2008;24:957-959
15. Pesando PM, Ghiringhello MP, Di Meglio G, Fanton G. Posterior chamber phakic intraocular lens (ICL) for hyperopia: Ten-year follow-up. *J Cataract Refract Surg.* 2007;33:1579-1584
16. Chung TY, Park SC, Lee MO, Ahn K, Chung ES. Changes in Iridocorneal Angle Structure and Trabecular Pigmentation with STAAR Implantable Collamer Lens During 2 Years. *J Refract Surg.* 2009;25:251-258
17. Allan B, Argeles-Sabate I, Mamalis N. Endophthalmitis Rates After Implantation of the Intraocular Collamer Lens: Survey of Users Between 1998 and 2006. *J Cataract Surg* 2009;35:766-769
18. Muzyka-Woźniak M, Ogar A. Anterior chamber depth and iris and lens position before and after phacoemulsification in eyes with a short or long axial length. *J Cataract Refract Surg.* 2016;42(4):563-568.



































STAAR Surgical AG Spain
Avenida Meridiana 216,
Despacho 4
Barcelona 08027
Spain

CE
0344

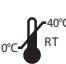

STAAR Surgical AG
Hauptstrasse 104
CH-2560 Nidau
Switzerland
Tel: +41 32 332 8888

SYMBOL GLOSSARY

| | | | |
|---|--|---|---|
|  | Medical device |  | Authorized representative in the European Community |
|  | Do not re-use |  | CE conformity marking per European Council Directive 93/42/EEC or European Council Regulation (EU) 2017/745 |
|  | Do not re-sterilize |  | Manufacturer |
|  | Do not use if the product sterile barrier system or its packaging is compromised |  | Date of manufacture |
|  | Body diameter (Optic diameter) |  | Country of manufacture – United States |
|  | Overall diameter |  | Country of manufacture – Switzerland |
|  | Single sterile barrier system with protective packaging outside |  | Unique Device Identifier |
|  | Use-by date |  | Catalogue number |
|  | Diopter |  | Right eye |
|  | Date |  | Left eye |
|  | Caution |  | Serial number |
|  | Contains biological material of animal origin |  | Spherical power |
|  | U.S. (Federal) law restricts this device to sale by or on the order of a physician |  | Cylindrical power |
|  | Store at room/ambient temperature. Do not freeze. Do not expose to temperature greater than 40°C |  | Axis |
|  | Health care center or Doctor |  | Spherical equivalent power |
|  | Sterilized using steam |  | Consult electronic instructions for use |

edfu.staar.com
+1-800-352-7842
+41 32 332 8888